



**BRAND  
NEW  
PRODUCTS!**



**OFFER  
PRICE  
£29.95**  
(RRP £44.25)

## Compoform UV

Compoform UV is a light curing composite for modelling.

- Can be combined with all types of waxes
- Burns out without leaving any residue
- Readily flowable yet stable consistency
- Uses standard UV lights

Can be used for:

- Modelling post and core restorations (extra oral)
- Modelling casting elements in the milling and attachment technique
- Extension of CoCr restorations
- Connecting bridges for soldering
- Blocking out undercuts




**OFFER  
PRICE  
£19.95**  
(each)

## Preparation burs & new shapes with SH cut

The new shape SH cut burs are suitable for quick removal of large quantities of plaster and composite. A smooth surface is achieved due to the relief cut technology.

- No chipping of plaster
- Chatter free processing
- Has proven relief cut technology
- Will last 3 times longer than normal burs
- The convex shapes allow quick removal of composite

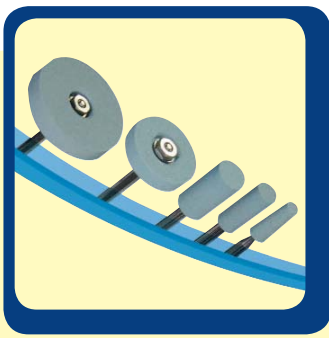

## Diagen

**Finishing Zirconia restorations has never been easier!**

Fast cutting, long lasting and the most aggressive sintered diamonds for zirconium and alumina oxide materials due to a special manufacturing process. Minimal pressure is all that is needed to achieve high grinding performance. 9 different shapes available!

**OFFER**  
buy any 3  
Diagens & get  
Fissure Tool  
Kit FREE!  
(RRP £41.95)



## Exacto-form

Dual component polyurethane resin for fabricating extremely precise dental models.

Due to its high edge stability, Exacto-form is perfectly suitable for precision-fit crown and bridge work.

Due to its very low viscosity, it is easily poured into impressions without entrapping air bubbles. This cost-effective model resin features absolutely precise reproduction of surfaces, dimensional stability and unequalled precision of fit with no expansion whatsoever.

Hardened resin can be drilled and trimmed with conventional burs.

Kits available in: Blue, grey, green, ivory and yellow

**OFFER**  
buy two kits  
of Exacto-form  
(£49.95 each)

receive **FREE** Micro-flex  
(34000141 worth £17.91)



## Pi-Ku-Plast HP36

For any task that requires a precision pattern resin.

**Fast 4 minutes set time, zero shrinkage and can be invested after 20 minutes!**

- Optimum contouring properties
- Available in 5 transparent colours
- Ideal for connecting bar sections for soldering
- Absorbs less monomer therefore less shrinkage
- Burnout is clean and complete
- During burnout it will not expand but collapses eliminating internal stress in the casting ring
- Available in 5 colours!

**OFFER**  
buy any kit  
(£99.95 each)

receive a **FREE** Preparation Bur  
(RRP £30.45)

